

Use of a High-resolution Melting Assay to Evaluate HIV Gag Region Diversity in HIV-infected Adults with Different Stages of HIV Infection



William I Towler¹, L Khaki¹, L Wang², M James¹, P Musoke^{3,4}, T Coates⁵, R Moore¹, G Kelen¹, O Laeyendecker^{1,6}, and SH Eshleman¹

¹Johns Hopkins Univ Sch of Med, Baltimore, MD, USA; ²Fred Hutchinson Cancer Res Ctr, Seattle, WA, USA; ³Makerere Univ-Johns Hopkins Univ Res Collaboration, Kampala, Uganda; ⁴Makerere Univ Sch of Med, Kampala, Uganda; ⁵Univ California at Los Angeles, Los Angeles, CA, USA; ⁶Lab of Immunoregulation, NIAID, NIH, Baltimore, MD, USA.

720 Rutland Ave.
Ross Bldg, Rm. 646
Baltimore, MD 21205
seshlem@jhmi.edu
Phone: 410-614-4734
Fax: 410-502-9244

Introduction

HIV diversity is usually studied by analyzing sequences from individual HIV variants. We developed a novel High Resolution Melting (HRM) assay to analyze HIV diversity without sequencing.

Methods

Template Preparation

PCR products used as templates in the HRM assay were generated using the ViroSeq HIV-1 Genotyping System (Celera, Alameda, CA) according to the manufacturer's instructions. ViroSeq PCR products were further diluted 1:10 to yield a concentration of ≈5 ng for use in the HRM assay.

HRM Assay

A region of HIV *gag* was amplified using the LightScanner Master Mix amplification buffer (Idaho Technology, Salt Lake City, UT). The resulting amplicons (~150-235 base pairs) were then analyzed on a LightScanner instrument (Model HR 96). The resulting data were processed to produce a melting curve for each sample representing the negative derivative of the data (-d[fluorescence]/d[temperature]). The left and right margins of each melting curve were marked at the positions at which the curve reached an angle of 30° relative to a horizontal baseline. The HRM score is defined as the difference (number of degrees) between the two borders.

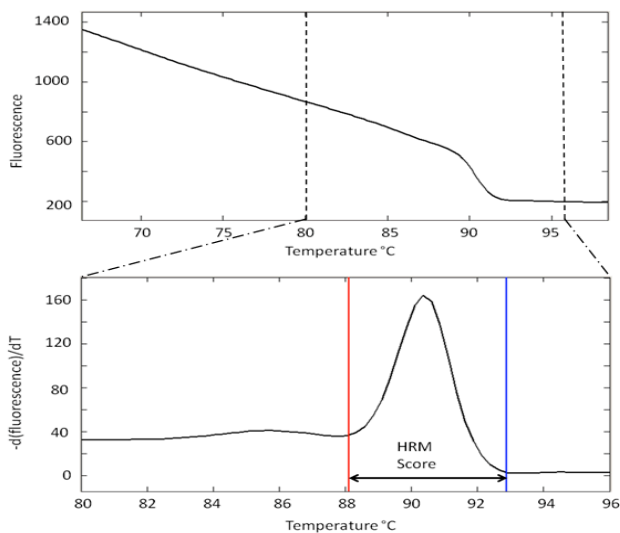


Figure 1. Analysis of HRM data to produce the HRM score.

Results

I. Proof of Concept

We first used the HRM assay to analyze plasmid controls and plasma samples from HIV-infected women and infants (9 mother-infant pairs). The median HRM score for infants (4.3, range 4.2-5.3) was higher than that for control plasmids (3.4, range 3.2-3.8, $p < 0.001$), and lower than that for mothers (5.7, range 4.4-7.7, $p = 0.005$, exact Wilcoxon rank sum test). The melting curves for each mother-infant pair are shown below (Figure 2).

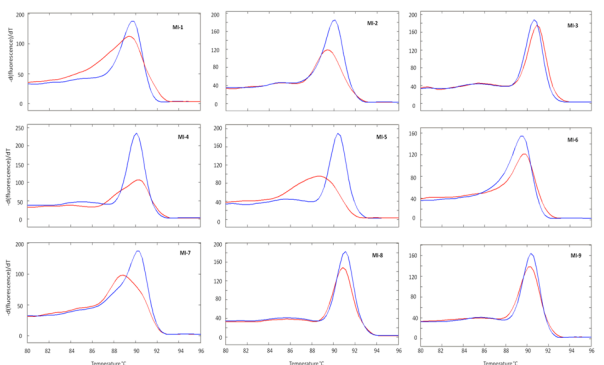


Figure 2. HRM data from HIV-infected women (red) and their infants (blue).

II. Assay Reproducibility

The same samples were analyzed with the HRM assay four times over a one year period. The intra-class correlation coefficient was 94% (95% CI: 89%-98%), indicating a high level of assay reproducibility.

III. Comparison of the HRM Assay to Traditional Sequence-based Diversity Measures

We compared HRM scores to traditional sequence-based measures of genetic diversity. For these studies, the HRM products from the 9 mother-infant pairs were cloned and sequenced (50 clones/woman, 20 clones/infant). Sequences were analyzed using MegaAlign, ClustalX2 and MEGA. Higher HRM scores were associated with higher genetic diversity ($p < 0.001$), higher complexity ($p < 0.009$), and higher Shannon entropy ($p < 0.022$). HRM scores were not associated with length variation ($p < 0.111$).

IV. Biologic Relevance

We compared the HRM scores of adults with different stages of HIV disease (Figure 3).
 • Acute (RNA positive, antibody negative, EXPLORE Study, $n = 20$)
 • Recent (median 187 days after last HIV negative test, EXPLORE Study, $n = 106$)
 • Chronic (HIV pos > 2 years, CD4 cell count > 50, JHU Moore Clinic and ED, $n = 38$)
 • AIDS (CD4 cell count < 50, JHU Moore Clinic and JHU ED, $n = 34$)

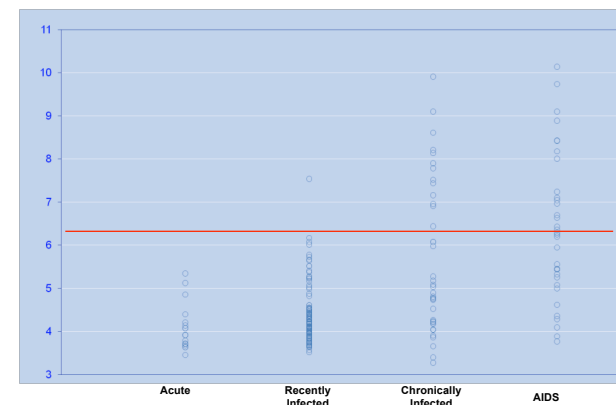


Figure 3. HRM scores from adults with different stages of HIV disease. The red line indicates an HRM score of 6.3 (mean + 3 SD of HRM scores for recently infected adults).

The median HRM score was lower in adults with acute HIV infection than in adults with recent HIV infection ($p = 0.0091$), and was lower in adults with recent HIV infection than in adults with chronic HIV infection and AIDS ($p < 0.0001$). The median HRM score was also lower in adults with chronic HIV infection than in adults with AIDS, but the difference was not statistically significant ($p = 0.061$). An HRM score of > 6.3 (mean + 3 SD for samples from adults with recent HIV infection) was associated with chronic HIV infection or AIDS ($p < 0.0001$, Fisher's exact test).

Conclusions

The HRM assay provides a novel, rapid method for assessing HIV diversity without sequencing. In this study, HIV *gag* region diversity increased as HIV infection progressed from acute HIV infection, to recent HIV infection, to chronic HIV infection and AIDS. This assay could be applied to any region of the HIV genome or to other genetic systems that exhibit DNA diversity.

Acknowledgements

The authors thank Jason McKinney, Justin Geiger, and Luke Stewart (Idaho Technologies) for helpful discussions.